


EXT10a/m/n

IMPLEMENTATION OF EXTERNAL ANALYSIS MODULES IN NOTOCORD-HEM™

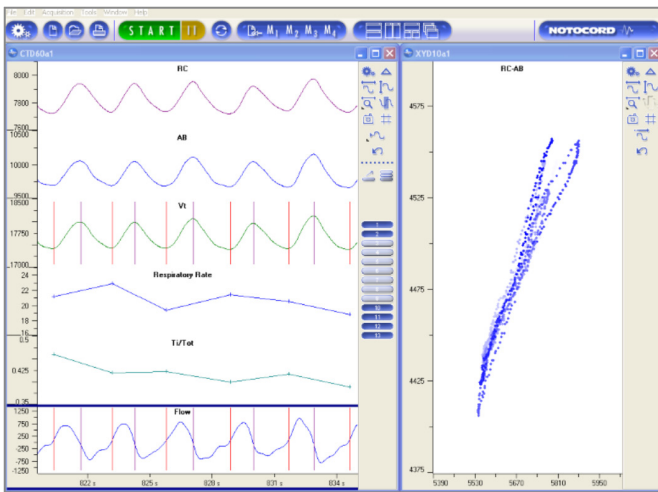

Respiration analyzer developed under EXT10a

KEY FEATURES

- Compatibility with Win32, Matlab® or .NET
- Integration of externally developed modules on NOTOCORD-hem™ platform
- Easy-to-create configuration window for external analysis modules
- Offline / online use of integrated modules
- Implementation by the user with or without NOTOCORD® involvement

CREATE YOUR OWN ANALYSIS MODULE

Customized analysis tools increase specificity in research studies such as detection of drug potential side effects.

Design and integrate in NOTOCORD-hem™ your own module developed under **Win32, MATLAB®** or **.Net!**

We created dedicated functionalities **EXT10a/m/n** that allow to implement external modules in our software platform to get as close as possible to the ideal analysis tool.

CONFIDENTIALITY AND CODE PROPERTY

Analysis modules based on specific modeling or innovative experiments remain confidential as they are encapsulated in a Dynamic Link Library (DLL) accessible to the owner only.

The source code developed for external analysis modules remains the property of the owner.

Our Product and Design teams can assist you in specifying **your requirements** and provide you with the appropriate information to develop your own analysis modules.


COMPATIBILITY WITH MATHWORKS MATLAB®

Processing and analysis of NOCOTORD-hem™ and **MATLAB®** data streams are possible in both environments using our library of functions **MPI10a**.

MPI10a allows retrieval of data streams from NOTOCORD-hem™ into MATLAB®, data reading and data writing from MATLAB® into NOTOCORD-hem™.

Original data streams are preserved. Raw data cannot be modified as all data written from MATLAB® are retrieved in separate module(s) in NOTOCORD-hem™.

NOTOCORD-hem™ Evolution is compatible with Mathworks MATLAB® from version R2006b and above.

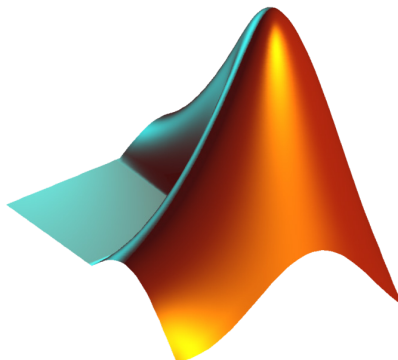

EXT10a/m/n

EUROPE - Headquarters

NOTOCORD SAS
Tel. +33 (0)1 34 80 00 00
information@notocord.com

NORTH AMERICA

NOTOCORD Inc.
Tel. +1 888 204 7770
information@notocord.com

JAPAN

PRIMETECH Corp.
Tel. +81 (0)3 3816 0851
sales@primetech.co.jp

SOUTH KOREA

IWOO Scientific Corp.
Tel. + 82 (0)2 3473 2332
sale@iwoo.co.kr

SEJONG-BIO
Tel. + 82 (0)4 2826 9473
sales@sejong-bio.co.kr

TAIWAN

MESONICS Systems Ltd.
Tel. + 886 2 2736 4066
info@mesonics.com.tw

CHINA

B&E Teksystems Ltd
Tel. + 86 10 8758 5773
bandetek@hotmail.com

INDIA

BIO-DOT Scientific, New Delhi
Tel: +91 11 2636 3490
upkar@biodotscientific.com

www.notocord.com
information@notocord.com
+33(0) 1 34 80 00 00

Working in a GLP environment?

21 CFR Part 11 compliance is achieved via our AccessManager security application, activating a highly secure and compliant Audit Trail.



Version 1.0 - © NOTOCORD - All rights reserved.

Because we have a policy of continuous product improvement, NOTOCORD reserves the right to change specifications without notice. The information provided in this brochure are not intended to form part of any offer or contract.

